

**3M<sup>™</sup> Imprint<sup>™</sup> 3** VPS Impression Material

# Accuracy at every step.

### **3M<sup>™</sup> Imprint<sup>™</sup> 3** VPS Impression Material

# No compromise on precision.

Imprint 3 VPS Impression Material's impressive balance of clinically relevant properties help you get a highly accurate impression on the first take. It is a versatile VPS material that performs at every step of your everyday procedures.



### Your procedure, your choice.

No matter what technique or case, you will find an Imprint 3 VPS Impression Material designed to suit the situation. Choose from a wide range of viscosities, from putty to light body, and two setting times.




### Catch every detail through superior hydrophilicity and flow behavior.

In moist conditions, an impression material's hydrophilicity and flowability are crucial to make void-free impressions and capture all details of the prepared teeth and margins.

Imprint 3 Wash Materials' superior hydrophilicity in the unset stage ensures great detail reproduction – even when blood and saliva are present.

To ensure you get the desired flow and optimal wetting of the tooth preparation, wash materials are available in three viscosities and two setting speeds for each.

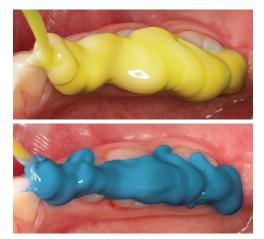

Flow to tooth and gingiva of 3M<sup>™</sup> Imprint<sup>™</sup> 3 Light Body (yellow) and 3M<sup>™</sup> Imprint<sup>™</sup> 3 Ultra-Regular Body (dark teal) wash materials. The more flowable 3M<sup>™</sup> Imprint<sup>™</sup> 3 Light Body creates a very smooth surface, whereas 3M<sup>™</sup> Imprint<sup>™</sup> 3 Ultra-Regular Body exhibits more structure.



# **Distortion-free removal** from the mouth due to overall **toughness**.

To avoid tearing upon mouth removal, a material needs both high tensile strength and high elongation potential. Toughness describes a combination of these two parameters.

Imprint 3 Wash Materials have outstanding tensile strength, elongation potential and overall toughness. Together, these features virtually eliminate the risk of impression distortion upon mouth removal.



# Fewer retakes thanks to recovery from deformation.

A material's memory is a function of its recovery from deformation and elongation. High recovery from elongation is especially important around the preparation margin.

Imprint 3 Material almost completely regains its original shape after elongation – reducing the risk of costly retakes.



### Your procedure, your choice.

3M<sup>™</sup> Imprint<sup>™</sup> 3 VPS Impression Material: Working and Setting Times

	Dispensing	Viscosity	Working Time at 23°C/74°F	Intra-oral Syringe Time	Intra-oral Setting Time
Product	System	low high	min:sec	min:sec	min:sec
Tray Materials					
3M <sup>™</sup> Imprint <sup>™</sup> 3 Penta <sup>™</sup> Putty			1:30		3:00
3M™ Imprint™ 3 Penta™ Heavy Body			2:00		3:30
3M <sup>™</sup> Imprint <sup>™</sup> 3 Penta <sup>™</sup> Quick Step Heavy Body			1:15		2:30
3M™ Imprint™ 3 Heavy Body			2:00		4:00
3M™ Imprint™ 3 Quick Step Heavy Body			1:15		2:30
3M™ Imprint™ 3 Monophase (Medium Body)			1:00		4:00
Wash Materials					
3M™ Imprint™ 3 Ultra-Regular Body			2:00	1:00	3:30*
3M™ Imprint™ 3 Quick Step Ultra-Regular Body			1:30	0:40	2:30
3M™ Imprint™ 3 Regular Body			2:00	1:00	3:30*
3M™ Imprint™ 3 Quick Step Regular Body			1:30	0:40	2:30
3M™ Imprint™ 3 Light Body			2:00	1:00	3:30*
3M™ Imprint™ 3 Quick Step Light Body			1:30	0:40	2:30

\*Intra-oral setting time is 3:00 min. when used in combination with 3M<sup>™</sup> Imprint<sup>™</sup> 3 Penta<sup>™</sup> Putty.

#### **Ordering Information**

#### 3M<sup>™</sup> Imprint<sup>™</sup> 3 VPS Impression Material

Refills — 3M <sup>™</sup> Penta <sup>™</sup> Materials	A		
Each refill includes: 2x 300ml base paste; 2x 60ml catalyst paste	Mixing U	es for 3M <sup>™</sup> Pentamix <sup>™</sup> 3 Automatic Init	
10771PH 3M <sup>™</sup> Imprint <sup>™</sup> 3 Penta <sup>™</sup> Heavy Body Refill   10971PH 3M <sup>™</sup> Imprint <sup>™</sup> 3 Penta <sup>™</sup> Quick Step Heavy Body Refill	P3808	3M <sup>™</sup> Imprint <sup>™</sup> 3 Penta <sup>™</sup> Heavy Body Cartridge	
10771PP 3M <sup>™</sup> Imprint <sup>™</sup> 3 Penta <sup>™</sup> Putty Refill Refills — Tray Materials	P3807	3M <sup>™</sup> Imprint <sup>™</sup> 3 Penta <sup>™</sup> Quick Step Heavy Body Cartridge	
Each refill includes: 4×50ml cartridges, 10 VPS mixing tips, green	P3806	3M <sup>™</sup> Imprint <sup>™</sup> 3 Penta <sup>™</sup> Putty Cartridge	
10771 3M <sup>™</sup> Imprint <sup>™</sup> 3 Heavy Body Refill	Accessories		
109713M <sup>™</sup> Imprint <sup>™</sup> 3 Quick Step Heavy Body Refill107813M <sup>™</sup> Imprint <sup>™</sup> 3 Monophase Medium Body Refill	71505	20x 3M <sup>™</sup> Intra-oral Syringes Green- Refill Pack	
Refills — Wash Materials Each refill includes: 4×50ml cartridges, 10 VPS mixing tips, yellow	71506	50x 3M <sup>™</sup> Intra-oral Syringes Green- Value Pack	
36956 3M <sup>™</sup> Imprint <sup>™</sup> 3 Ultra-Regular Body Refill	77580	1x 3M <sup>™</sup> Garant <sup>™</sup> Cartridge Dispenser 1:1/2:1	
36955 3M <sup>™</sup> Imprint <sup>™</sup> 3 Quick Step Ultra-Regular Body Refill	71512	50x 3M <sup>™</sup> Penta <sup>™</sup> Mixing Tips- Red	
10779 3M <sup>™</sup> Imprint <sup>™</sup> 3 Regular Body Refill	71450	50x VPS mixing tips, green	
10979 3M <sup>™</sup> Imprint <sup>™</sup> 3 Quick Step Regular Body Refill	71452	50x VPS mixing tips, yellow	
10778 3M <sup>™</sup> Imprint <sup>™</sup> 3 Light Body Refill	71460	50x VPS intra-oral tips, bite, clear	
10978 3M <sup>™</sup> Imprint <sup>™</sup> 3 Quick Step Light Body Refill	71462	50x VPS intra-oral tips, yellow	
	7307	1x VPS Tray Adhesive 17ml	



3M Oral Care 2510 Conway Avenue St. Paul, MN 55144-1000 USA

Phone 1-800-634-2249 Web 3M.com/dental

#### Warranty

3M warrants this product will be free from defects in material and manufacture at the time of purchase. 3M MAKES NO OTHER WARRANTIES INCLUDING ANY IMPLIED WARRANTY OF MERCHANTABILITY OR FITNESS FOR A PARTICULAR PURPOSE. User is responsible for determining the suitability of the product for user's application. If this product is defective within the warranty period, your exclusive remedy and 3M's sole obligation shall be repair or replacement of the 3M product.

#### Limitation of Liability

Except where prohibited by law, 3M will not be liable for any loss or damage arising from this product, whether direct, indirect, special, incidental or consequential, regardless of the theory asserted, including warranty, contract, negligence or strict liability.

3M, Garant, Imprint, Penta and Pentamix are trademarks of 3M or 3M Deutschland GmbH. Please recycle. © 3M 2021. All rights reserved.