



1.25 Gallon

BioSonic[®] UC125

Ultrasonic Cleaning System

Relax, it's a BioSonic.

New LCD technology offers customized operation

Largest tank capacity in its class

Operates at a very quiet 63 dB


coltène 
whaledent[®]

BioSonic® UC125

Largest Tank Volume in its Class



Multi-tasking

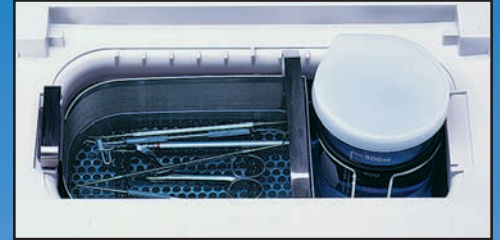

Full basket



Multi-beaker



Useful Features


Faster Draining

The UC125 drains in approximately 45 seconds (half the time of the UC100XD)!


C-Power

An exclusive feature that ensures constant power against varying solution levels and instrument loads.

Tank Size (Gallons)

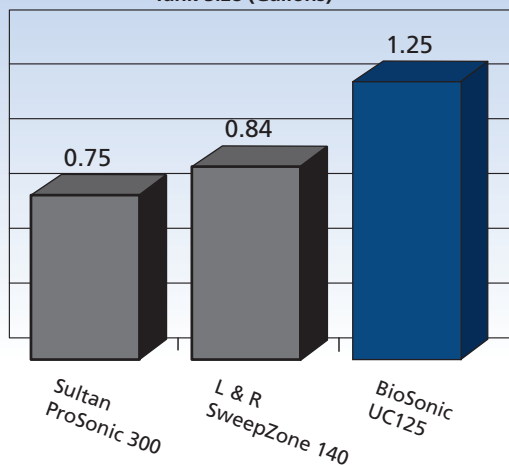


60% Larger than other "In Class" cleaners

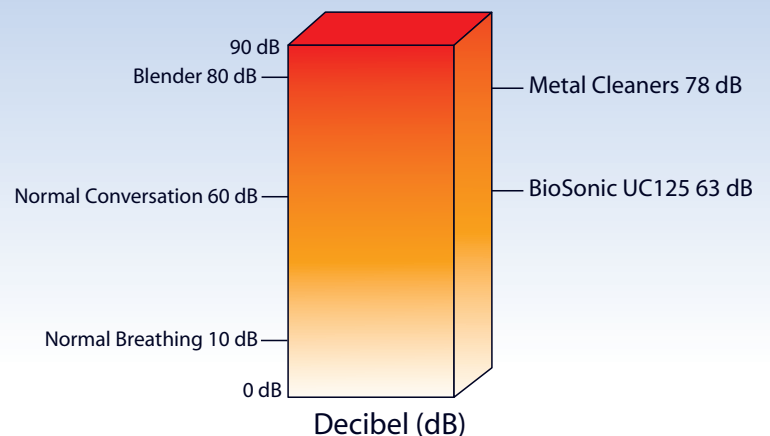

20% quieter than the competitors' metal ultrasonic cleaners


Now with **LCD Technology** for customized operation

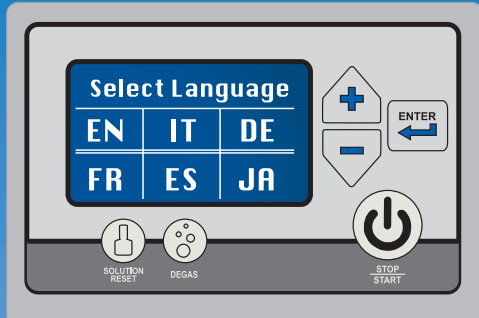


MULTIPLE LANGUAGES
Choose one of six different languages



DEGAS BUTTON
Easily prepare your solution for maximum effective cleaning

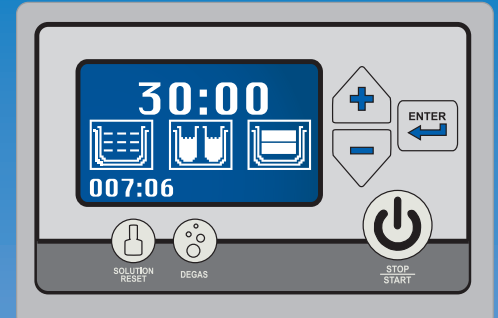

SELECTABLE RUN TIMES
Have the freedom to customize cleaning time by operation

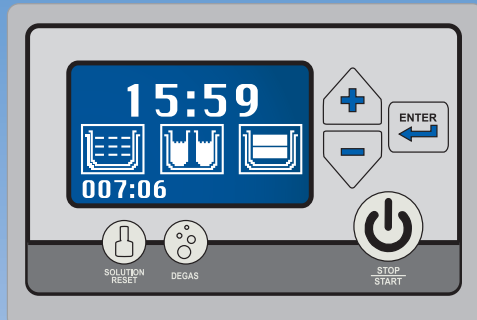

COUNTDOWN TIMER
Know exactly when your instruments will be ready for the autoclave



SOLUTION TRACKING
Stop guessing how long your solution has been in use

UC125 Accessories


Cross-flow baskets offer enhanced instrument cleaning, 25% deeper with easy lift-out handles.

UC151XD



Full-Size Basket

UC152XD

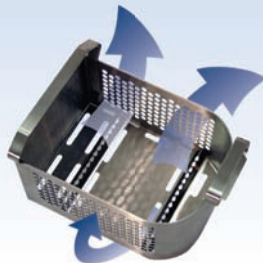

3/4 Size Basket

Hu-Friedy IMS®*



Cassette Rack

Hu-Friedy IMS®*



Instrument Cassette

* IMS is a registered trademark of Hu-Friedy Inc.

BioSonic® UC125 Ultrasonic Cleaning System

Technical Data

Tank Capacity:	1.25 gallons (4.8 liters)
Use Capacity:	1.00 gallons (3.8 liters)
Overall Dimensions:	16.5 x 10.5 x 13.75 inches (42 x 26.7 x 35 cm)
Tank Dimensions:	11.25 x 5.62 x 5 inches (28.6 x 14.3 x 12.7 cm)
Output:	137 watts
Operating Frequency:	52 kHz

For maximum performance and longevity of your BioSonic Ultrasonic Cleaner, we recommend the use of genuine BioSonic solutions and accessories.

Order Information

UC125	BioSonic Ultrasonic Cleaner Accessories not included
UC151XD	Full Size FingerGuard™ Instrument Basket
UC152XD	Three-Quarter FingerGuard Instrument Basket
UC155XD	Accessory Kit Contains (1) UC152XD, (2) UC153, and (2) UC53
UC53	Borosilicate Beaker 600 ml With cover and positioning ring
UC153	One Beaker Holder
UC164XD	Cassette Rack for Hu-Friedy® Space Saver Seven Holds 2 cassettes

BioSonic® Solution Concentrates in MeterDose® bottles



Solution:	UC30 BioSonic Super Concentrate General Purpose Cleaning Solution	UC32 BioSonic Enzymatic Ultrasonic Cleaning Concentrate	UC34 BioSonic Plaster and Stone Remover Powder Concentrate	UC38/UC42 BioSonic Germicidal Ultrasonic Cleaning Concentrate	UC39 BioSonic Ultrasonic Cement Remover Concentrate
Contents:	16 oz. (473 ml)	8 oz. (236 ml)	1.85 lb. (840 g)	8 oz. (236 ml) UC38 16 oz. (473 ml) UC42	16 oz. (473 ml)
Diluted Yield:	8 gal. (30 liters)	8 gal. (30 liters)	3 gal. (11.4 liters)	8 gal. (30 liters) UC38 16 gal. (60 liters) UC42	.4 gal. (1.5 liters)
Dilution Ratio:	1:64	1:128	1:135	1:128	1:2
Contaminants:	Blood, debris, saliva and general cleaning	Protein soils (dried blood, tissue, saliva) and general cleaning. Effectively breaks down and loosens hardened debris.	Plaster and stone	Blood, debris, saliva and general cleaning. Cleans instruments as it disinfects the ultrasonic bath.	Temporary cements and permanent cements (including zinc-phosphate, polycarboxylate)
Items:	Instruments, crowns, bridges, dentures and ultrasonic scaler tips	Instruments, crowns, bridges, dentures and ultrasonic scaler tips	Dentures, partials, castings, denture flasks, mixing bowls and non-aluminum articulators	Instruments, crowns, bridges (no plastic parts)	Crowns, veneers, brackets and instruments

Coltène/Whaledent Inc.
235 Ascot Parkway
Cuyahoga Falls, OH 44223 / USA
Tel. USA & Canada 1 800 221 3046
+1 330 916 8800
Fax +1 330 916 7077
info@coltenewhaledent.com

Coltène/Whaledent AG
Feldwiesenstrasse 20
9450 Altstätten / Switzerland
Tel. +41 (0)71 757 53 00
Fax +41 (0)71 757 53 01
office@coltenewhaledent.ch

Authorized EU Representative:
Coltène/Whaledent GmbH + Co. KG
Raiffeissenstrasse 30
89129 Langenau / Germany
Tel. +49 (0)7345 805 0
Fax +49 (0)7345 805 201
productinfo@coltenewhaledent.de

Coltène/Whaledent website:
<http://www.coltenewhaledent.com>

coltène
whaledent®