

Clinpro™ Prophylaxis Paste and Angles.



Clinpro™

Prophylaxis Paste

- 3M ESPE's Clinpro™ prophylaxis paste removed stains as well as NUPRO® Coarse, while being significantly less abrasive to both dentin and enamel.
- Using Clinpro prophylaxis paste results in less gloss reduction which means the high esthetics and luster of microfill restorations will continue to enhance your patient's smile.

Prophylaxis Angles

- Specifically designed to meet the needs of hygienists, the single-use disposable Clinpro prophylaxis angles are available with soft or firm non-latex* cups.

* This product is made of a synthetic material and does not intentionally contain natural rubber latex.

"Clinpro Prophylaxis Paste is an outstanding alternative for cleaning and polishing both natural teeth, direct, and indirect restorations."

Trish Jones, RDH, BS, Technical Advisor for Aurum Ceramic © LVI

Give your patients the best of both worlds — a great polish with less abrasion.

Why choose between a coarse or fine grit prophylaxis paste when you can use both? The unique perlite particles within Clinpro™ prophylaxis paste are formulated to break down as you polish from a coarse to fine grit. This unique feature provides the stain removal of a coarse paste with less abrasion to the enamel and dentin.

New




Your patients receive an excellent polish and will be delighted with the result.

Clinpro prophylaxis paste also incorporates features found in many other prophylaxis pastes:

- Contains fluoride
- Neutral pH value
- Refreshing mint and bubble gum flavors
- Available in unit dose cups


After 15 Seconds

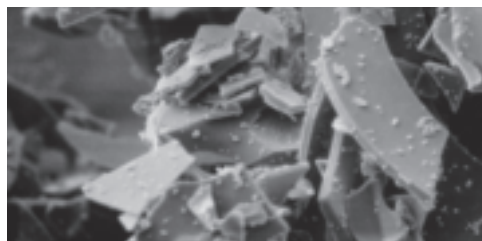

Large cleaning particles within the paste (magnified 3000x) break into fine polishing particles in a matter of seconds.

The benefits show. Compare ClinproTM prophylactic paste to the leading pastes.

Quality Performance from 3M ESPE

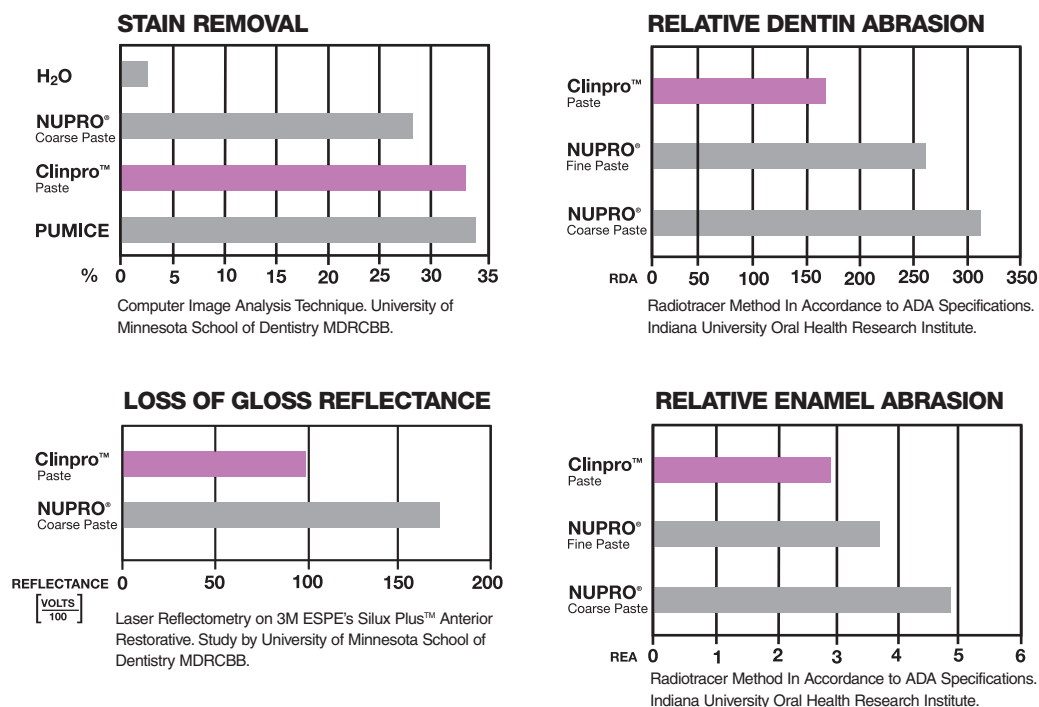

Ordering Information



Item Number	Product Information
12611	ClinproTM Prophylactic Paste, Mint Flavor Includes 200 unit dose cups and one prophylactic ring.
12612	ClinproTM Prophylactic Paste, Bubble Gum Flavor Includes 200 unit dose cups and one prophylactic ring.
12600	ClinproTM Prophylactic Angles, Firm Cup Includes 144 single-use disposable angles with firm non-latex* cups
12601	ClinproTM Prophylactic Angles, Soft Cup Includes 144 single-use disposable angles with soft non-latex* cups
12626	ClinproTM Sealant Introductory Kit, Syringes Contains: 2 – 1.2ml Clinpro Sealant Syringes; 1 – 3ml 3M TM ESPE TM Scotchbond TM Etchant Syringe; 20 – Sealant Syringe Tips (Black); 25 – Etchant Syringe Tips (Blue)
12621	ClinproTM Sealant Introductory Kit, Bottles Contains: 2 – 6ml Clinpro Sealant Bottles; 1 – 9ml 3M TM ESPE TM Scotchbond TM Etchant Bottle; 2 – Brush Handles; 60 – Brush Tips; 50 – Fiber Tips; 1 – Mixing Well with Cover

For a free sample pack of Clinpro Prophylactic Paste/Angles (5 cups mint flavor paste, 5 cups bubble gum flavor paste and 3 angles), call 1-888-640-7121 ext. 1.

*This product is made of a synthetic material and does not intentionally contain natural rubber latex.

Technical Hotline: 1-800-634-2249 www.3MESPE.com

3M ESPE

Dental Products

3M Center
Building 275-2SE-03
St. Paul, MN 55144-1000
USA

3M Canada

Post Office Box 5757
London, Ontario N6A 4T1
Canada
1-800-265-1840 ext. 6229



10% post-consumer waste paper

Printed in USA
© 3M 2004

3M, ESPE, Clinpro, Scotchbond, and Silux Plus are trademarks of 3M or 3M ESPE AG. NUPRO is a registered trademark of Dentsply International. Used under license in Canada.

70-2009-3705-3