

**3M** Science.  
Applied to Life.™

**3M™ Imprint™ 3**  
VPS Impression Material

# Accuracy at every step.

Imprint 3 VPS Material's impressive balance of clinically relevant properties helps you capture a highly accurate impression on the first take. It is a versatile VPS material that performs at every step of your everyday procedures.



## Captures every detail.

The wash materials' superior hydrophilicity and flow allow for void-free impressions with great reproduction — even when blood and saliva are present.



## Distortion-free mouth removal.

Excellent overall toughness helps Imprint 3 impressions resist tearing upon mouth removal.



## Good memory for fewer retakes.

Imprint 3 material's memory allows impressions to regain their original shape after stretching during mouth removal. This almost 100% recovery reduces the chance of costly remakes.





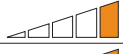

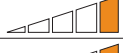

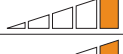







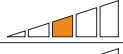








## 3M™ Imprint™ 3 VPS Impression Material

- Available in a broad range of viscosities and setting times
- For 1-step and 2-step impression techniques



# Your procedure, your choice.

## 3M™ Imprint™ 3 VPS Impression Material: Working and Setting Times

Product	Dispensing System	Viscosity low  high	Working Time at 23°C/74°F min:sec	Intra-oral Syringe Time min:sec	Intra-oral Setting Time min:sec
<b>Tray Materials</b>					
3M™ Imprint™ 3 Penta™ Putty			1:30		3:00
3M™ Imprint™ 3 Penta™ Heavy Body			2:00		3:30
3M™ Imprint™ 3 Penta™ Quick Step Heavy Body			1:15		2:30
3M™ Imprint™ 3 Heavy Body			2:00		4:00
3M™ Imprint™ 3 Quick Step Heavy Body			1:15		2:30
3M™ Imprint™ 3 Monophase (Medium Body)			1:00		4:00
<b>Wash Materials</b>					
3M™ Imprint™ 3 Ultra-Regular Body			2:00	1:00	3:30*
3M™ Imprint™ 3 Quick Step Ultra-Regular Body			1:30	0:40	2:30
3M™ Imprint™ 3 Regular Body			2:00	1:00	3:30*
3M™ Imprint™ 3 Quick Step Regular Body			1:30	0:40	2:30
3M™ Imprint™ 3 Light Body			2:00	1:00	3:30*
3M™ Imprint™ 3 Quick Step Light Body			1:30	0:40	2:30

\*Intra-oral setting time is 3:00 min. when used in combination with 3M™ Imprint™ 3 Penta™ Putty.

## Ordering Information

### 3M™ Imprint™ 3 VPS Impression Material

Item no.	Product information	Item no.	Product information
<b>Refills — 3M™ Penta™ Materials</b>		<b>Cartridges for 3M™ Pentamix™ 3 Automatic Mixing Unit</b>	
Each refill includes: 2x 300ml base paste; 2x 60ml catalyst paste			
10771PH	3M™ Imprint™ 3 Penta™ Heavy Body Refill	P3808	3M™ Imprint™ 3 Penta™ Heavy Body Cartridge
10971PH	3M™ Imprint™ 3 Penta™ Quick Step Heavy Body Refill	P3807	3M™ Imprint™ 3 Penta™ Quick Step Heavy Body Cartridge
10771PP	3M™ Imprint™ 3 Penta™ Putty Refill	P3806	3M™ Imprint™ 3 Penta™ Putty Cartridge
<b>Refills — Tray Materials</b>		<b>Accessories</b>	
Each refill includes: 4x50ml cartridges, 10 VPS mixing tips, green			
10771	3M™ Imprint™ 3 Heavy Body Refill	71505	20x 3M™ Intra-oral Syringes Green-Refill Pack
10971	3M™ Imprint™ 3 Quick Step Heavy Body Refill	71506	50x 3M™ Intra-oral Syringes Green-Value Pack
10781	3M™ Imprint™ 3 Monophase Medium Body Refill	77580	1x 3M™ Garant™ Cartridge Dispenser 1:1/2:1
<b>Refills — Wash Materials</b>		71512	50x 3M™ Penta™ Mixing Tips- Red
Each refill includes: 4x50ml cartridges, 10 VPS mixing tips, yellow		71450	50x VPS mixing tips, green
36956	3M™ Imprint™ 3 Ultra-Regular Body Refill	71452	50x VPS mixing tips, yellow
36955	3M™ Imprint™ 3 Quick Step Ultra-Regular Body Refill	71460	50x VPS intra-oral tips, bite, clear
10779	3M™ Imprint™ 3 Regular Body Refill	71462	50x VPS intra-oral tips, yellow
10979	3M™ Imprint™ 3 Quick Step Regular Body Refill	7307	1x VPS Tray Adhesive 17ml
10778	3M™ Imprint™ 3 Light Body Refill		
10978	3M™ Imprint™ 3 Quick Step Light Body Refill		

#### Warranty

3M warrants this product will be free from defects in material and manufacture at the time of purchase. 3M MAKES NO OTHER WARRANTIES INCLUDING ANY IMPLIED WARRANTY OF MERCHANTABILITY OR FITNESS FOR A PARTICULAR PURPOSE. User is responsible for determining the suitability of the product for user's application. If this product is defective within the warranty period, your exclusive remedy and 3M's sole obligation shall be repair or replacement of the 3M product.

#### Limitation of Liability

Except where prohibited by law, 3M will not be liable for any loss or damage arising from this product, whether direct, indirect, special, incidental or consequential, regardless of the theory asserted, including warranty, contract, negligence or strict liability.

3M, Garant, Imprint, Penta and Pentamix are trademarks of 3M or 3M Deutschland GmbH. Please recycle.  
© 3M 2021. All rights reserved.



3M Oral Care  
2510 Conway Avenue  
St. Paul, MN 55144-1000 USA  
Phone 1-800-634-2249  
Web 3M.com/dental