

ACTIVA™ BioACTIVE-CEMENT™

PULPDENT®

Inspired by Nature



The BioACTIVE CONNECTION - Natural Lute, Bond & Seal

ACTIVA is Different

Bioactivity stimulates mineral apatite crystal formation that penetrates and fills micro-gaps, knits together the restoration and the tooth, guards against recurrent caries, and seals margins against microleakage and failure – nature's way.

APATITE FORMATION SEALS AND PROTECTS TEETH

ACTIVA stimulates the natural remineralization process. Apatite crystals form ionic bonds that connect the restoration to the tooth, fill micro-gaps, and seal margins at the restorative-tooth interface, defining a new standard of care.

ABSORBS SHOCK AND STRESS

ACTIVA's patented rubberized-resin provides greater resistance to chipping and fracture than any other dental cement.

CAD/CAM ADVANTAGE

ACTIVA BioACTIVE-CEMENT forms strong bonds to zirconia, ceramic, lithium disilicate, and metal. Ideal for CAD/CAM and implant dentistry. Tissue compatible.

FAST AND EASY TO USE

No etching. No bonding. Dual cure. Tack cure for easy removal of excess.

DURABLE AND INSOLUBLE


ACTIVA's bioactive resin matrix does not wash out or crumble. It is moisture tolerant and ideal for the moist oral environment.

RESPONDS TO pH CHANGES

ACTIVA is a dynamic material and the only one with continuous release and recharge of calcium, phosphate and fluoride for long-term benefits.

NO BISPHENOL A

ACTIVA is safe for all your patients. It contains no Bisphenol A, no Bis-GMA and no BPA derivatives.

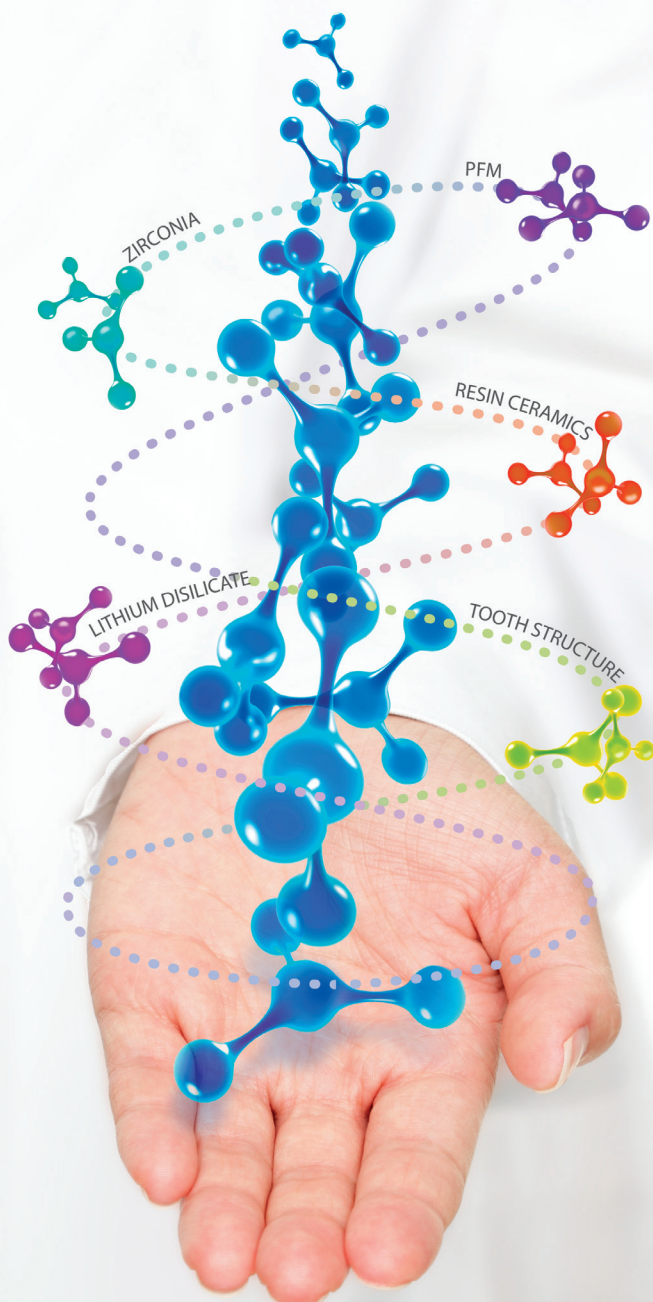

ACTIVA™ BioACTIVE-CEMENT™

SINGLE PACK: 5mL/7gm syringe + 20 automix tips
(15 straight black + 5 with bendable metal cannula)

VC1A2	A2 Opaque	\$114.99
VC1T	Translucent	\$114.99

VALUE PACK: 2 x 5mL/7gm syringes + 40 automix tips
(30 straight black + 10 with bendable metal cannula)

VC2A2	A2 Opaque	\$206.99
VC2T	Translucent	\$206.99

Mineral Apatite Formation SEM Analysis of ACTIVA™ BioACTIVE-CEMENT™ Surface after 21 Days in Saline

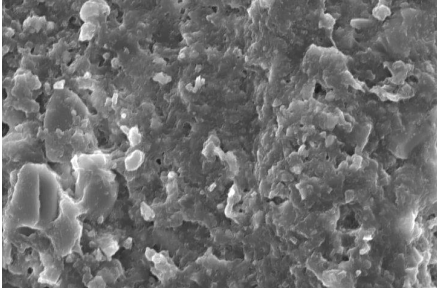


Fig 1 ACTIVA BioACTIVE-CEMENT Control, no saline 3000x

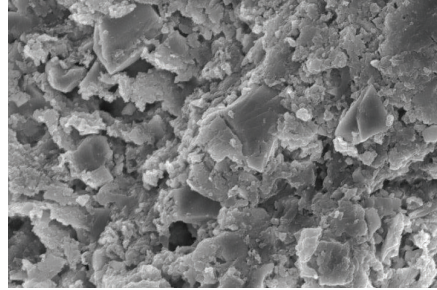



Fig 2 ACTIVA BioACTIVE-CEMENT 21 days in saline 3000x

Compared to the no saline control, scanning electron microscope (SEM) imaging and energy-dispersive X-ray spectroscopy (EDS) after 21 days in saline shows significant increase in calcium and phosphorus ion concentrations, and decrease in carbon and silica ions, indicating that mineral apatite deposits are forming on the surface.

A simple cementation procedure: self-adhesive, syringe delivery, no trituration, easy clean up



1 Tooth is prepared to receive a crown. Note retentive crown prep.

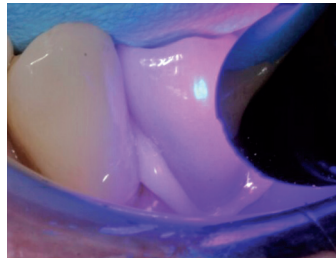

2 Crown filled with ACTIVA BioACTIVE-CEMENT is seated and tack cured 1-2 seconds.






3 Excess cement is easily removed



4 Shows finished case

Photos courtesy of Dr. Frankie Shull

Physical Properties of ACTIVA™ BioACTIVE-CEMENT™ Compare Favorably with Leading Cements

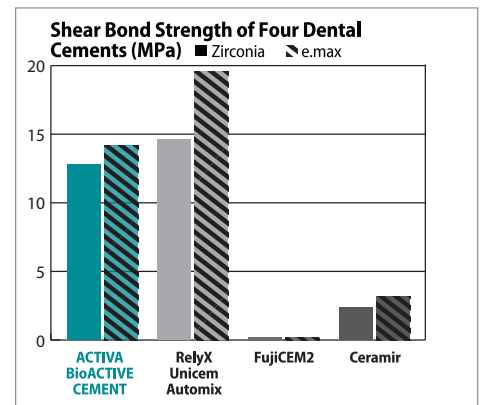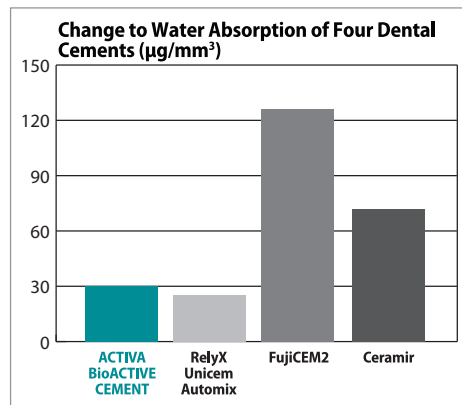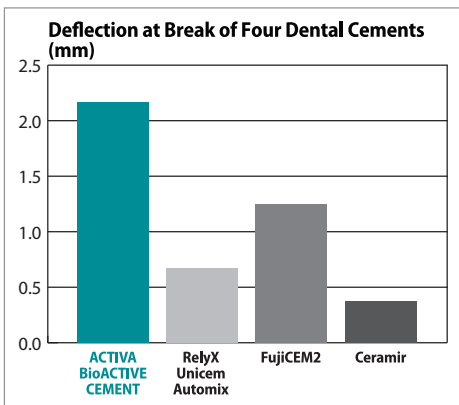

RelyX, FujiCEM and Ceramir are trademarks of 3M ESPE, GC and Doxa respectively

Contact PULPDENT or your local PULPDENT sales representative about ACTIVA today!